


Neo-Vent™

Infant T-piece Resuscitation Circuit

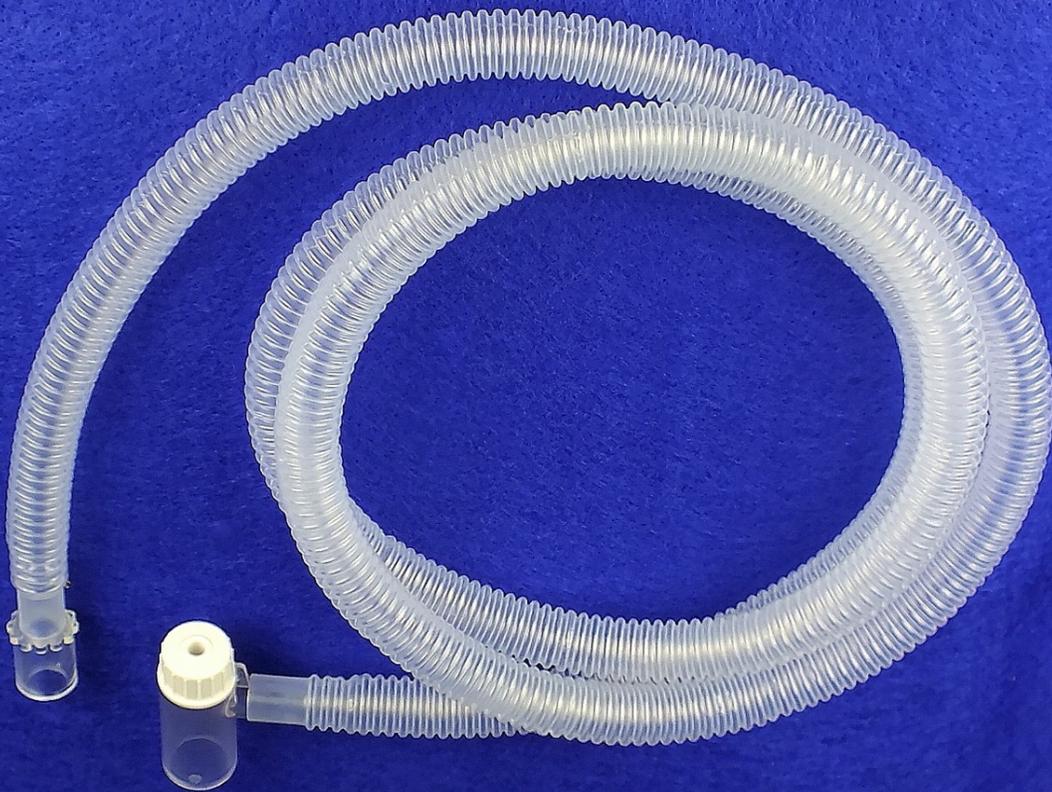

Westmed
When Your Care Is Critical

Neo-Vent™

Infant T-piece Resuscitation Circuit



When newborn resuscitation efforts are undertaken, product performance should never be in doubt. The Neo-Vent™ infant T-piece resuscitation circuit is a simple, dependable product that is ready when you need it.

Features & Benefits

- ✓ Compatible circuit for use on Fisher & Paykel Healthcare Neopuff™ 900-series Infant T-piece Resuscitator
- ✓ For use on infants weighing <10 Kg
- ✓ Universal T-piece connector: 15 mm ID
- ✓ Circuit hose length: 60 inches
- ✓ Variable PIP adjustable from 0 to 45 cm H₂O
- ✓ Operator adjusts PEEP with knob on T-piece and controls respiratory rate and inspiratory time with thumb pressure on T-piece valve
- ✓ Less hand-fatigue; no required compression of a resuscitator bag
- ✓ Single-patient Use, disposable

The Neo-Vent infant T-piece circuit complies with standard practice guidelines recognized by ILCOR and AHA's Neonatal Resuscitation Programs.



0417 Gas Supply Tubing Assembly

Part Number	Description	Case Qty
0415	Infant T-piece Resuscitation Circuit	10
0416	Infant T-piece Resuscitation Circuit & Gas Supply Tubing	10
0417	Gas Supply Tubing Assembly	10
0420	Infant T-piece Resuscitation Circuit with Infant Mask	10
0421	Infant T-piece Resuscitation Circuit with Neonatal Mask	10



Distributed by:

Westmed, Inc.
5580 S. Nogales Highway
Tucson, AZ 85706
Toll Free: (800) 975-7987
Fax: (520) 294-6061
www.westmedinc.com



Visit www.westmedinc.com
for a complete listing of
Westmed products